



# iM70

## Veterinary Monitor

- 12.1-inch high-resolution color TFT-LCD screen
- Full touch screen
- Customizable shortcuts for simple and intuitive operation
- No-fan design reduces dust and creates a no noise environment
- Defibrillation and electrosurgical interference protection
- Basic parameters monitoring
- Optional capnography and anesthesia monitoring available
- Suitable for feline, canine and other animals

Standard Parameters: 3/5-lead ECG, HR, RESP, EDAN SpO2, EDAN NIBP, PR, 2-TEMP

Optional Parameters: Nellcor OxiMax™ SpO2, SunTech NIBP, 2-IBP, EDAN G2 CO2 (Sidestream), Respiration CO2 (Mainstream and Sidestream), Masimo Anesthesia Gas/O2 (iM70 VET)

**1200**

NIBP  
Measurement

**120 h**


Trend  
Review

**60**

Alarm  
Review

**120 s**

Frozen  
Waveform

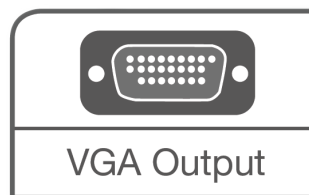



iM70

12.1" Color TFT-LCD Touch Screen



iM8


10.4" Color TFT-LCD Touch Screen




No-fan Design



Pitch Tone

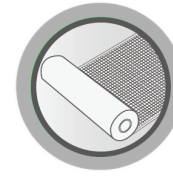

Thermal Recorder

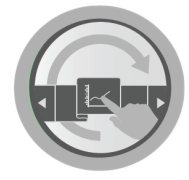

Shortcut Menu

## End Tidal Carbon Dioxide Monitoring



### EDAN G2 CO<sub>2</sub> (sidestream)

Superior water trap design for accurate monitoring

iCARBTM algorithm with intelligent CO<sub>2</sub> pseudo wave identification technology

Multiple sampling accessories as options

### Respironics CO<sub>2</sub> (mainstream/sidestream)

Plug & play module design

Dehumidification tube instead of water trap

Low sampling rate of 50ml/min

## MASIMO Anesthetic Gas/O<sub>2</sub> (Mainstream/Sidestream)



Unique plug & play module for mainstream anesthetic gas monitoring

The no moisture sampling line design suitable for all clinical application

A stable low flow of 50ml/min

## Multiple Display Modes



Standard Display



Large Font



Trend View



OxyCRG



Back to Bed



Wall Mount



Rolling Stand