



iM70 Veterinary Monitor

- 12.1-inch high-resolution color TFT-LCD screen
- Optional Full touch screen
- Customizable shortcuts for simple and intuitive operation
- No-fan design reduces dust and creates a no noise environment
- Defibrillation and electrosurgical interference protection
- Basic parameters monitoring
- Optional capnography and anesthesia monitoring available
- Suitable for feline, canine and other animals

Standard Parameters: 3/5-lead ECG, HR, RESP, EDAN SpO2, EDAN NIBP, PR, 2-TEMP

Optional Parameters: Nellcor OxiMax™ SpO2, SunTech NIBP, 2-IBP, EDAN G2 CO2 (Sidestream), Respiration CO2 (Mainstream and Sidestream), Masimo Anesthesia Gas/O2 (iM70 VET)

1200

NIBP
Measurement

120 h

Trend
Review

60

Alarm
Review

120 s

Frozen
Waveform





iM70
12.1" Color TFT-LCD Screen




iM60
10.4" Color TFT-LCD Screen




No-fan Design



Pitch Tone

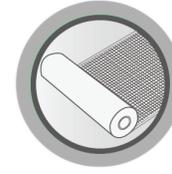

Thermal Recorder
(Optional)

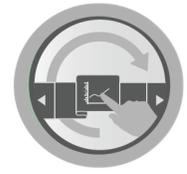

Shortcut Menu

End Tidal Carbon Dioxide Monitoring



EDAN G2 CO2 (sidestream)

Superior water trap design for accurate monitoring
iCARBTM algorithm with intelligent CO2 pseudo wave identification technology
Multiple sampling accessories as options

Respironics CO2 (mainstream/sidestream)

Plug & play module design
Dehumidification tube instead of water trap
Low sampling rate of 50ml/min

MASIMO Anesthetic Gas/O2 (Mainstream/Sidestream)



Unique plug & play module for mainstream anesthetic gas monitoring
The no moisture sampling line design suitable for all clinical application
A stable low flow of 50ml/min

Multiple Display Modes



Standard Display



Large Font



Trend View



OxyCRG



Back to Bed



Wall Mount



Rolling Stand